
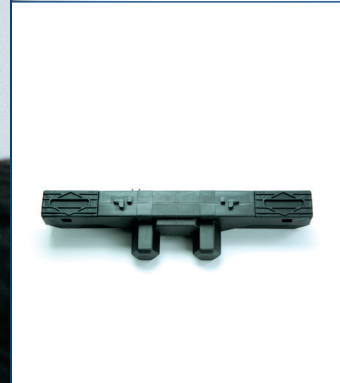


INVACARE®  
**VALUE CARE**  
SEMI-ELECTRIC BED

---

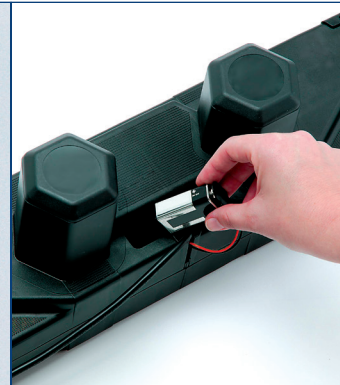

---

The **Invacare® Value Care Semi-Electric Bed** utilizes a single motor design, making the bed lightweight and easy to assemble. It features the popular universal bed ends with the durable U.S.-made Invacare bed frame. The advanced pendant design allows the head and knee sections to move either simultaneously or independently at the push of a button.

Invacare beds, mattresses and rails, excluding therapeutic support surfaces and bariatric beds, when used as a system and sold new and unused after August 1, 2007, meet the voluntary FDA Guidelines for bed rail entrapment. Visit **[www.invacare.com/bedsafety](http://www.invacare.com/bedsafety)** for additional safety and entrapment information.

---

## Features

U.S.-made bed deck is the industry-leading design with riveted joints to provide years of maintenance-free service

Single motor design makes the bed lightweight and easy to assemble

9V battery backup provides convenient protection in case of power outages



Advanced pendant design allows head and knee sections to operate simultaneously or independently

Universal bed ends are durable and impact and scratch resistant

Semi-electric operation:

- Electric head and knee sections
- Manual hi/low operation

## HCPCS Code: E0295

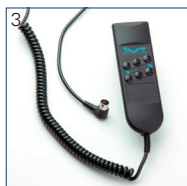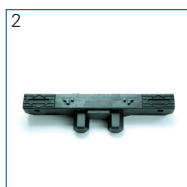

1. 9V battery
2. Single motor
3. Hand pendant



**Invacare Corporation**  
www.invacare.com

### USA

One Invacare Way  
Elyria, Ohio  
44035-4190  
(800) 333-6900

### Canada

570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

©2015  
Invacare Corporation.  
All rights reserved.  
Trademarks are  
identified by the  
symbols ™ and ®.  
All trademarks are  
owned by or licensed  
to Invacare Corporation  
unless otherwise  
noted. Specifications  
are subject to change  
without notification.  
Form No. 06-024  
Rev. 08/15  
1506106



Model #	Dimensions	Sleep Surface	Weight Capacity	Shipping Weight	Carton Dimensions	Regulation Compliance	Limited Warranty
VC5310	88" L 15"-23" H 36" W	80" L 36" W	Safe Working Load: 450 lb.* Max. Patient Weight: 350 lb.	N/A	N/A	UL962 CSA Certified	5 Years Welds 2 Years Mechanical/ Electrical
VC5000 Headspring				45 lb.	43" L 8" W 36" D		
VC5890 Footspring				81 lb.	47" L 9" W 36" D		
5301IVC Universal Bed End				55 lb.	39" L 7" W 24" D		

\*Includes patient, mattress and all accessories.

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.