

INVACARE® FULL-ELECTRIC BED



The **Invacare® Full-Electric Bed** offers the greatest convenience for the patient and caregiver. The easy-to-use pendant control provides motorized positioning of the upper body and knees, and also adjusts bed frame height for efficient use of accessories and for safe and easy transfer.

Invacare beds, mattresses and rails, excluding therapeutic support surfaces and bariatric beds, when used as a system and sold new after August 1, 2007, meet the voluntary FDA Guidelines for bed rail entrapment. Visit www.invacare.com/bedsafety for additional safety and entrapment information.

Features

High impact bed end panels are more durable, impact and scratch resistant; washable for easy cleaning

Universal bed ends can be used as either the headboard or footboard on the bed

Rubber caster tread for better tracking with non-marring material

DC low voltage, compact motor system
Quiet motor operation and fast cycle time on lowering height of bed for efficient bed set-up

Motor system with waterproof enclosures is rated IP66. Washable with standard water pressure. The bed must not be power washed

Color-coded connections on junction box for easier maintenance


Double-insulated junction box provides extra protection against electrical shock

Ergonomic hand pendant comes with durable strain relief and heavy-duty cord. Waterproof casing for easy cleaning, conveniently attaches to bed sheet with clip or rail with hand pendant cradle (sold separately)

Hand pendant cradle available for convenient storage of hand pendant model # 26975

Interchangeability simplifies storage, reduces inventory complexity

HCPCS Code: E0297

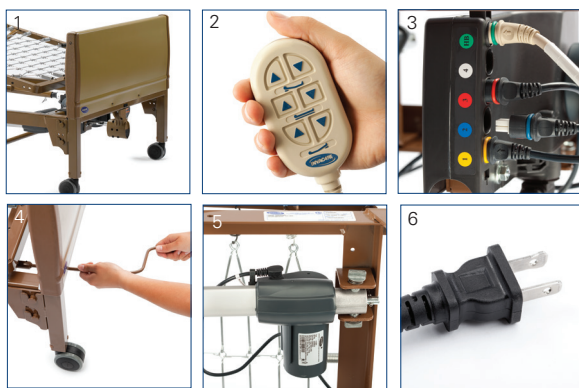

1. Patented universal bed ends
2. Ergonomic hand pendant
3. Color coded connections
4. Emergency manual crank
5. Thermal overload protection
6. Double insulated electrical components/wiring



Invacare Corporation
www.invacare.com

USA

One Invacare Way
Elyria, Ohio
44035-4190
(800) 333-6900

Canada

570 Matheson Blvd. E.,
Unit 8
Mississauga, Ontario
L4Z 4G4, Canada
(800) 668-5324

©2015

Invacare Corporation.
All rights reserved.
Trademarks are identified by the symbols™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.
Form No. 10-123
Rev. 08/15
1506101



Model #	Dimensions	Sleep Surface	Weight Capacity	Shipping Weight	Carton Dimensions	Regulation Compliance	Limited Warranty
5410IVC	88" L 15"-23" H 36" W	80" L 36" W	Safe Working Load: 450 lb.* Max. Patient Weight: 350 lb.	N/A	N/A	UL962 CSA Certified	5 Years Welds 2 Years Mechanical/ Electrical
5000IVC Headspring				45 lb.	43" L 8" W 36" D		5 Years Welds 2 Years Mechanical/ Electrical
5490IVC Footspring				81 lb.	47" L 9" W 36" D		5 Years Welds 2 Years Mechanical/ Electrical
5301IVC Universal Bed End				55 lb.	39" L 7" W 24" D		

*Includes patient, mattress and all accessories.

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.